


AUTOMATIC DONOR COUCH



AUTOMATIC DONOR COUCH

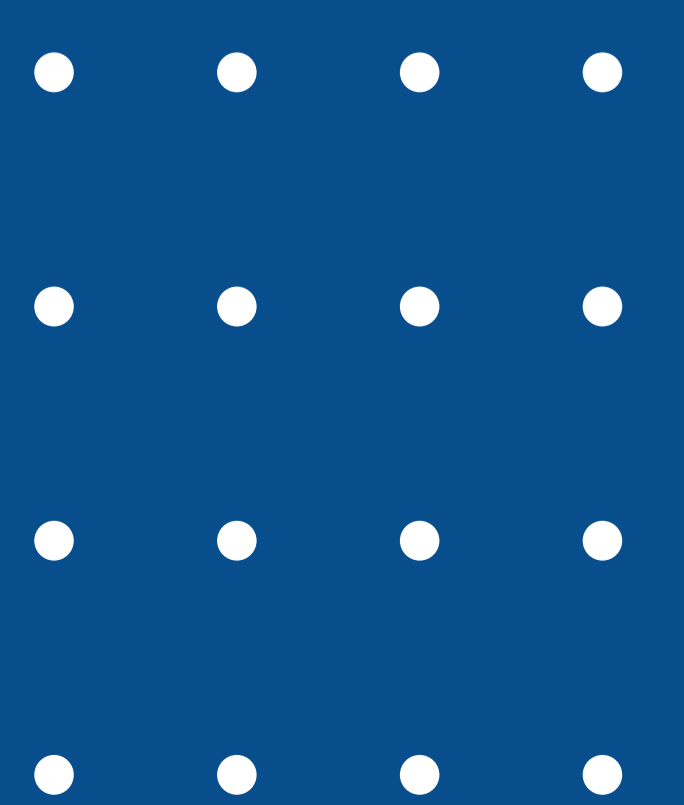

ELEKTRO Dual actuator based blood donor couch is specially designed to make blood withdrawal easier, safe and provide comfortable position of the donor. Donor's position can be easily adjusted to form either a comfortable easy chair or a bed. All surfaces are clean and smooth, easily accessible and can be kept clean. Electrical and mechanical accessories are housed inside the structure making the product safer. Seat elements are designed with ergonomic consideration to maximize donor comfort.

FEATURES

- Two Geared Actuator to adjust head, back rest and foot end to Ensure safety and comfort to the donor. Heavy Duty Tilt adjustment can be done by using remote control, as well as manual with battery support (optional) in case of power failure.
- Specially designed to Provide a comfortable position for the donor.
- Variable position of head rest, back rest, foot rest and arm rest can be achieved through actuator to change donor's posture. (Donor's head can be lowered immediately and legs can be lifted above heart level so that blood can flow back to the brain and other vital organs, in case of vasovagal attack)
- All frame made of heavy gauge mild steel square pipe with duly powder coated.
- 4 Lockable castors for easy mobility
- Comfortable chair type with 4 inch soft cushioning mattress washable upholstery and comfortably wide arm-rests.

TECHNICAL SPECIFICATION

Model No.	EDC-02
Power supply	230 VAC
Weight	60kg
Lifting Capacity	150kg
Control	Corded Remote
Power Consumption	100 W
Movement Actuation	DC Motor
Size in Normal Condition	62 x 205 x 144 cm x (W x D x H)
Size In Stretched Condition	62 x 205 x 110 cm x (W x D x H)
Length of Arm Rest	60cm
Width of Arm Rest	15cm
Upholstery	Soft Upholstery of 4.0 thick
Battery Support	-----